


Turn a Duo Coronet to Mono Settings

1. Connect the coronet to the App


Open the Well Red app.

Select *Coronet*.


Select the **cogwheel** just under the green bar – the settings icon.

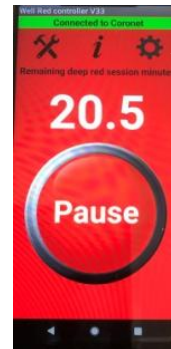

Select *Therapist zone*

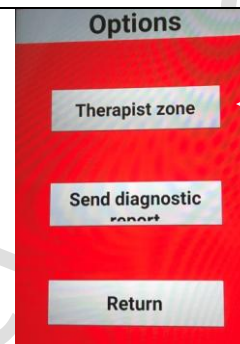

The screen asks for a code.

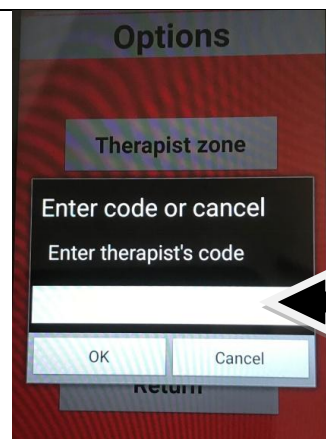


To change from Duo to Mono settings, enter:

`#Update,S1440,s1`

Select OK.

A few seconds later, your Coronet will blink once.

This means that the new Mono setting is in place.

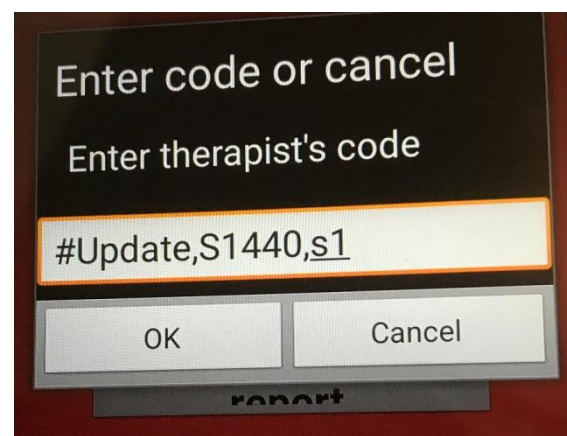

Turn a Duo Coronet to Mono Settings


If you want to reset your Duo from the Mono settings back to the Duo, you can enter this code:

`#Update,S720,s720`

Select OK.

A few seconds after you selected OK, your Coronet will blink once.

This means that the Duo setting is back in place.

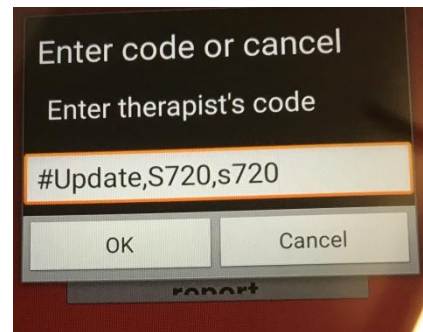

Well Red Pty Ltd.